



Dry Acid Concentrate Mix System

Fast. Reliable. Compact. User-Friendly.

Rockwell Medical's dry acid concentrate mixer is fast, compact, and easy-to-use. Our mixer is engineered for dialysis facilities seeking a streamlined, cost-effective solution for onsite mixing.

510(k) approved for exclusive use with Rockwell Medical's CitraPure® and Dri-Sate® dry acid powders, our mixer simplifies operations while upholding high quality standards and reliability.

Quick Installation, User-Centric Design

1 Rapid Mixing

Mixes 100 gallons in approx. 3 minutes¹ – over 90% faster than traditional systems.

2 48" Angled Side Port


Ergonomic design for easy, mess-free powder loading, cleaning and maintenance.

3 Compact Footprint

Fits through standard doors. Ideal for space-conscious facilities.

4 Plug-and-Play Installation

Installs quickly. Simply connects to your clinic's existing power and RO water for consistent, high-quality performance.

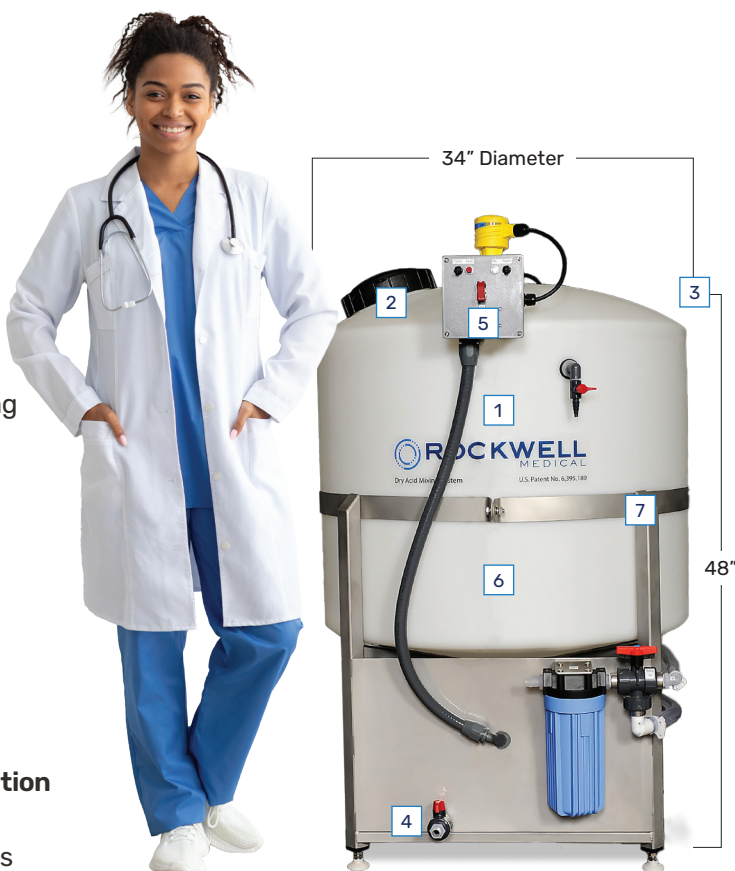

5 One-Touch Operation

Simple interface. Just press the button – mixer stops automatically when ready to use.

6 No Onsite Technician Required

Self-installation and easy maintenance. No manufacturer service contracts required.

7 Made in the U.S.A.

Durable, corrosion-free stainless steel construction for long-term reliability. Fast turnaround from order to delivery.

¹Data on File. Available Upon Request. Source: Rockwell Medical, Inc.

Why Choose Rockwell Medical Over Other Mixers?

FEATURE:	ROCKWELL MEDICAL'S MIXER	OTHER MIXERS
Mix Time	~3 minutes	>45 minutes
Installation	Simple, standard plumbing, various voltages	Complex, bulky setup
Operation	One-touch, walk-away ready	Multi-step with prompts
Maintenance	No service tech required, non-proprietary parts	Requires manufacturer support and service tech
Build	Marine-grade (316) stainless steel	Proprietary, complex components
Size	Compact design, small footprint 48" x 34"	Large footprint

Implementation & Ongoing Support

- Mix full 100 gallons batch using 4 cases of Rockwell Medical dry acid (CitraPure® or Dri-Sate®).
- No onsite technician required.
- Conduct specific gravity testing using a hydrometer and graduated cylinder (ranges provided in installation guide).
- If initial batch passes, solution is ready to use.
- Ongoing mixer support provided by Rockwell Medical.

Recommended Accessories (available from Rockwell Medical)

ITEM #	DESCRIPTION
PS-311	Hydrometer (C Series), Ranges 1.180 – 1.250
PS-24774-750	500 ml Graduated Cylinder
PS-GX05-9	Hytrex Replacement Cartridge Filter 5 micron
D4-Series or DC-Series	CitraPure Dry Acid Concentrate (25-gal mix)
	Dri-Sate Dry Acid Concentrate (25-gal mix)

Experience the Rockwell Medical Difference

Exceptional Quality, Excellent Customer Care, Reliable Delivery

Rockwell Medical develops, manufactures, commercializes, and distributes a portfolio of hemodialysis products for dialysis providers worldwide. We manage our own fleet of delivery trucks through our subsidiary, Rockwell Transportation, to ensure timely delivery and personalized service. Rockwell Medical has global reach with a growing network of partners and distributors in more than 30 countries. We are the #1 reliable supplier of hemodialysis concentrates. When other suppliers can't fulfill their commitments, you can rely on Rockwell Medical to deliver on-schedule and on-budget, every time.

Choose Reliability. Choose Simplicity. Switch to Rockwell Medical.

Contact Sales Support:

800.449.3353
sales@rockwellmed.com
www.rockwellmed.com

